



# EASYCAP

EEG Recording Caps and Related Products

EASYCAP GmbH  
Steingrabenstrasse 14  
DE-82211 Herrsching  
Germany  
Delivery Address: Am Anger 5, DE-82237 Woerthsee-Etterschlag

Tel +49 (0) 8153 88702-00  
Fax +49 (0) 8153 88702-10  
[www.easycap.de](http://www.easycap.de)  
[info@easycap.de](mailto:info@easycap.de)


## 74Ch actiCAP Snap, layout "CAC-64"/ M1

Cap with holders for acticap slim electrodes, additional holders acc. M1 layout

- Green Holders: Label 1-32, Hardware Channel 1 – 32      ● Yellow Holders: Label 1-32, Hardware Channel 33 - 64  
● Blue holder: Label & hard-wired REF      ● Black holder: Label & hard-wired GND      ○ additional white holders

### Electrode Names and Number Labels

**actiCAP**  
CCCsnap

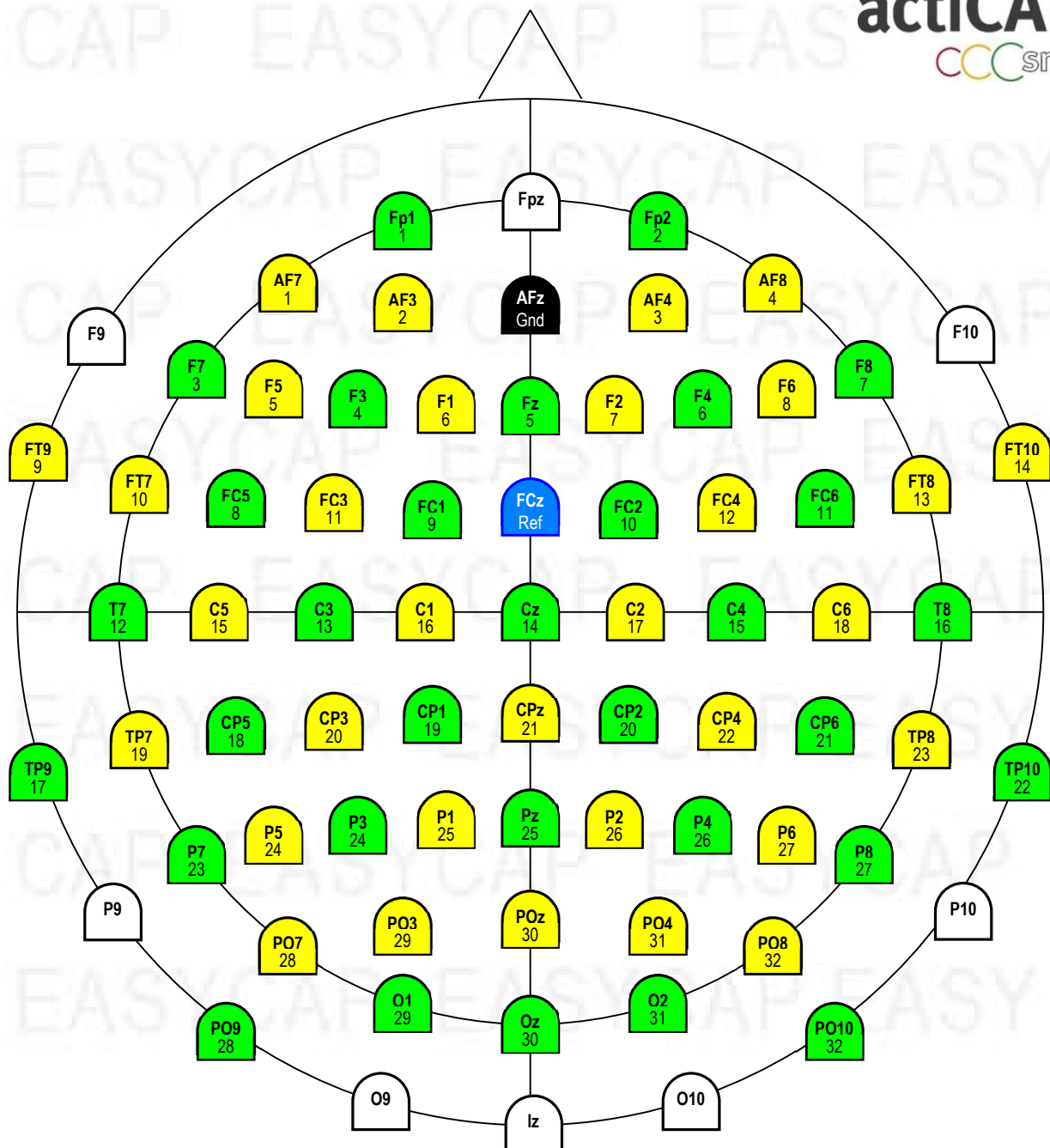

## Details



### Ordering Information

For ordering please give **Article Number, Cap Cut, and Size**

(Example: CACS-74-X1, Caucasian, 58):

- Article Number: **CACS-74-X1**
- Cap Cut: **Caucasian** or **Asian**
- Size (given in cm head circumference):
  - Adult caps: **52, 54, 56, 58, 60, 62, 64** (average male: 58, average female: 56)
  - Children caps: **50** (3-4 years), **52** (5-10 years), **54** (11-14 years)
  - Infant caps: **42, 44** (7 month), **46, 48** (2 years)

The catalogue-number comprises the cap as described, serial number, and this document; all packed in a labelled cardboard box. For further information about accessories or consumables please visit our website or contact our local distributor.

### Cap

White Subtemporal Cap with integrated chin belt.

All sizes made from High Comfort Fabric

Options: *Caucasian* or *Asian*, *Size*

The holders are mounted according to the figure. The positions are chosen from 10%-System. Positions TP9/TP10 approximately meet left and right mastoids. Additional holders in the white positions.

Apart from using the cap holders, it is also possible to use some electrodes on bare skin, e.g. for recording EOG. To do so, the electrodes can be attached to the skin with double-sided adhesive rings (washers). The grip of the washer is better, when both, skin and electrode, are degreased with alcohol prior to application.

### Electrode Assignment

The actiCAP electrodes come in bundles of 32, each labelled 1-32.

Plug the GND-Electrode into the black holder at AFz.

Plug the REF electrode into the blue holder at FCz.

Plug the electrodes 1-32 of the first bundle into the green holders (label 1-32), for amplifier channels 1-32.

Plug the electrodes 1-32 of the second bundle in the yellow holders (label 1-32), for amplifier channels 33-64.

The velcro straps attached to the cap help you to organize and fixate the cable wires.

## Supplementary Files

Each cap comes with a memory stick containing two bvef-files, one for use with REF, one for NO\_REF. Use these bvef-files to assign the electrode names to channel numbers while creating the BrainVision Recorder Workspace. The bvef-file also contains the coordinates of the electrodes. In addition, the stick contains a bvef-readme file, a guide map to illustrate the sequence of electrode numbers, and this document.

## Accessories

The cap can be fixed either by the integrated chin belt alone, or additionally / alternatively by a chest belt (not included). Normally the decision between these two methods can be made according to what is more comfortable for the subject. Under difficult conditions the chest belt will hold the cap in place a little better.

## Service

In case of any further questions

- Regarding consumables (e.g. electrolyte gel, syringes) and their usage, regarding cap, holders, layout, belts, please contact our local distributor or [info@easycap.de](mailto:info@easycap.de). Please have the serial number of your cap at hand (syntax CACSxxxxx, located on the backside label of the cap).
- Regarding actiCAP electrodes and associated electronics and software, please contact the local Brain Products Distributor or [techsup@brainproducts.com](mailto:techsup@brainproducts.com).